

APSW 4-Ply, Wire-Reinforced, Mandrel-Wrapped, Silicone Suction Hose Platinum-Cured; Smooth or Convoluted O.D.


ow volatile grade, platinum-cured silicone hose for critical pharmaceutical, biopharm, cosmetic, and food applications comprises the liner of AdvantaSil APSW. The liner is slipped onto a mandrel and wrapped with polyester mesh fabric, 302/304 stainless steel reinforcing wire and additional silicone. It is then cured into a homogeneous hose, enhancing its pressure and vacuum capabilities. AdvantaSil APSW has undergone extensive physical, chemical and biological testing and meets USP Class VI, FDA CFR 177.2600, ISO 10993 and European Pharmacopoeia 3.1.9 standards.




www.advantapure.com

KEY FEATURES



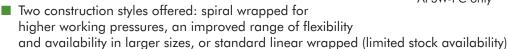

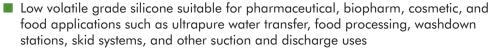

- Hardness value of 70 Shore A
- Rated for full vacuum
- Crush and kink resistant; durable and heavy-duty
- Linear wrapped style is listed by the National Sanitation Foundation for food equipment materials (NSF-51)
- Liner certified free of animal-derived ingredients
- General temperature range: -100°F (-73.3°C) to 350°F (176.6°C)
- Resists temperature extremes, chemical attack, ozone, radiation, moisture and environmental exposure
- Sterilizable by autoclave, CIP and SIP processes
- Documented lot traceable with identification on packaging
- Documented quality control per ISO 9001:2015

Complete validation package available on request








Liner

APSW-P & APSW-PC only



APSW 4-Ply, Wire-Reinforced, Mandrel-Wrapped Silicone Suction Hose

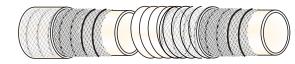

SPECIFICATIONS

儿		<u> </u>												
					Wo	rking	Bu	rst						
	Product	I.	D.	0	.D.	Vacuum		ssure	Pres	sure N	Iinimu n	n We	ight S	tandard
	Number	(in.)	(mm)		(mm)	Rating		0°F	at 7		Bend			Length
				(refe	rence)	(in./Hg.)		1°C)	,		Radius			(ft.)
							(PSI)	(Bar)	(PSI)	(Bar)	(in.)	(lb)	(kg)	
	APSW-SP	_ S1	1 00ТН	O D	Spin	AT W	D A DD	ΕD						
					/				1.500	1040	1.05	001	410	10
	APSW-SP-0500	.500	12.70	.910		>29	380	26.2		104.8		.281	.418	12
	APSW-SP-0750	.750	19.05	1.16	29.46	>29	315	21.7	1260	86.9	1.75	.359	.534	12
	APSW-SP-1000		25.40	1.41	35.81	>29	250	17.2	1000	68.9	3.0	.438	.652	12
	APSW-SP-1500		38.10	1.91	48.51	>29	250	17.2	1000	68.9	4.0	.685	1.019	12
	APSW-SP-2000		50.80	2.41	61.21 73.91	>29	205 185		820 740	56.5	6.0	.875	1.302	12
	APSW-SP-2500	2.500	76.20	2.91	86.61	>29	190	12.8		51.0	10.5 CF	1.050	1.830	12
	APSW-SP-3000 APSW-SP-4000	4.000	101.60	3.41	112.01	>29	140	13.1 9.7	760 560	52.4 38.6	CF		2.545	12
	AP3W-3P-4000	4.000	101.60	4.41	112.01	>29	140	9.7	360	30.0	Сг	1.710	2.545	12
APSW-SPC — CONVOLUTED O.D., SPIRAL WRAPPED														
	APSW-SPC-0500	.500	12.70	.950	24.13	>29	400	27.6	1600	110.3	1.0	.281	.418	12
	APSW-SPC-0750	.750	19.05	1.20	30.48	>29	338	23.3	1350	93.1	1.5	.359	.534	12
	APSW-SPC-1000	1.000	25.40	1.45	36.83	>29	288	19.9	1150	79.3	2.5	.438	.652	12
	APSW-SPC-1500	1.500	38.10	1.95	49.53	>29	250	17.2	1000	68.9	3.0	.685	1.019	12
	APSW-SPC-2000	2.000	50.80	2.45	62.23	>29	205	14.1	820	56.5	5.5	.875	1.302	12
	APSW-P -							ED —		E SUI		LAST	1	
	APSW-P-0750	.750	19.05	1.16		>29	150	10.3	600	41.4	2.5	.359	.534	13
	APSW-P-1000	1.000	25.40	1.41	35.81	>29	150	10.3	600	41.4	3.5	.438	.652	13
	APSW-PC	-C	ONVOL	0.D.	AR V	RAPI	PED -	- WH	LE SU	PPLIE	ES LAS	ST		
	APSW-PC-0500	.500	12.70	.950	24.13	>29	150	10.3	600	41.4	2.0	.208	.309	13
	APSW-PC-0750	.750	19.05	1.20	30.48	>29	150	10.3	600	41.4	2.5	.307	.457	13
	APSW-PC-1000	1.000	25.40	1.45	36.83	>29	150	10.3	600	41.4	3.5	.390	.580	13
	APSW-PC-1500	1.500	38.10	1.95	49.53	>29	150	10.3	600	41.4	4.0	.625	.930	13
	APSW-PC-2000	2.000	50.80	2.45	62.23	>29	125	8.6	500	34.5	6.0	.775	1.153	13

Nominal wall dimensions for all sizes of APSW-SP and APSW-P are .205" or 5.21mm; for APSW-SPC and APSW-PC: .225" or 5.72mm. Custom lengths and hose assemblies with an assortment of fittings are available. CF = Consult Factory

NOTE: When products are used as part of an assembly, the pressure ratings of fittings may be less than hose pressure ratings above. Please consult your sales representative.

Address: NewAge Industries AdvantaPure®

145 James Way

Southampton, PA 18966 U.S.A.

Toll Free Phone: 888-755-4370
Phone: 215-526-2151
Toll Free Fax: 888-258-4293
Fax: 215-526-2167

Web Site: www.advantapure.com E-mail: sales@advantapure.com Advanta Püre®



DISTRIBUTED BY:

Fluitek Componentes de Processos Ltda Rua Santo Bassan 111 - Bosque das Palmeiras CEP 13086-745 Campinas - SP Fone (19) 3257.3400 e-mail: atendimento@fluitek.com.br -

www.fluitek.com.br

